

veterinary technology for life

Eickemeyer®

www.eickemeyer.com

Veterinary Small Animal Dentistry Catalogue



Telephone +1 519 273 5558



4808146



4808258



4808259



4802037

Profident Mobile Dental Unit

The Profident Dental Unit is the result of years of experience offering applications for both general practitioners and dental specialists. Compact, mobile, adjustable height and clean design allow for easy sanitization before and after use and proper positioning for the operator.

The Profident offers the following features:

- High-speed handpiece with push button friction grip for burrs
- Low speed air motor and straight handpiece for polishing and heavy cutting
- Three-way syringe (air/water/mist)
- Foot control including water switch

Item No.	
4808146	Profident Dental Unit
4808258	NSK Low Speed Dental Handpiece <ul style="list-style-type: none"> • Up to 22,000 rpm • Air motor included • Internal water cooling
4808259	NSK Low Speed Dental Handpiece – Straight (HP) <ul style="list-style-type: none"> • 1 to 1 Direct drive • Up to 25,000 rpm • Single port spray
4802037	NSK Low Speed Dental Handpiece – Contra Angle (RA) <ul style="list-style-type: none"> • 1 to 1 Direct drive • Up to 40,000 rpm • Single port spray
4804743	Low Speed contra-angle <ul style="list-style-type: none"> • Fits air motor above
4808151	Extra 2 l water bottle

Eickemeyer®

More than 6,000
products online
www.eickemeyer.com



© Eric Issele – shutterstock.com

Profident PLUS Mobile Dental Unit

The Profident PLUS Dental Unit, in addition to all the features of the Profident, adds a built-in piezo dental scaler with LED illumination in the handpiece dramatically increasing your work efficiency.

The Profident PLUS adds offers the following:

- Built-in piezo dental scaler with LED illumination, variable pitch and water flow to more easily remove tartar and calculus
- All Profident features

Item No.

48008148	Profident PLUS Dental Unit
4808259	NSK Low Speed Dental Handpiece – Straight (HP) <ul style="list-style-type: none"> • 1 to 1 Direct drive • Up to 25,000 rpm • Single port spray
4802037	NSK Low Speed Dental Handpiece – Contra Angle (RA) <ul style="list-style-type: none"> • 1 to 1 Direct drive • Up to 40,000 rpm • Single port spray
4802627	Low Speed straight handpiece <ul style="list-style-type: none"> • Including air motor
4804289	Flat pointed scaler tip
4804290	Flat wide scaler tip
4804291	Flat pointed scaler tip
4800578	Flat pointed P1T scaler tip
4804743	Low Speed contra-angle <ul style="list-style-type: none"> • Fits air motor above
4808151	Extra 2 l water bottle



48008148



4808258



4808259



4802037



4804289



4804290



4804291



4800578



4808147

Dental Unit Air Compressor

Item No.

- 4808147 Dental Unit Air Compressor
- Required for Profident (4808146)/ Profident PLUS (4808148) Dental Unit
 - Mountable on the base of Profident/ Profident PLUS
 - Quiet, 40 dB (A) (as a fan)
 - Low maintenance
 - Compact design
 - Powerful: 340 watt
 - 230 V, 8 bar/116 psi, 37 l/min
 - 6 l capacity
 - Dimensions (in cm): L 45 x W 43 x D 20



4802062

High Speed Dental Handpiece with LED

Suitable for Profident/Profident PLUS Dental Units.

Item No.

- 4802062 High Speed Dental Handpiece with LED
- Interchangeable burr with push button
 - Rotation speed: 300,000 rpm
 - 32,000 lux LED light
 - No power supply required and fits most dental units



4809125



4802238



4809126



4809127



4809128

Small Animal Dental Burr Set

Item No.		Quantity	Form
4809122	Complete set, consisting of:		
4809125	Dental Burr (HP) • Size 23	2	round
4802238	Dental Burr (FG) • Size 23	5	round
4809126	Dental Burr (FG) • Size 10	5	round tapered
4809127	Dental Burr (FG) • Size 12	5	round tapered
4809128	Dental Burr (FG) • Size 16	5	round tapered

Dental Burr Set (RA) and Polishing Cups (RA)**Item No.**

- 4801447 Dental Burr Set (RA)
- Wide selection of burs to accommodate all your routine procedures on dogs and cats
 - Set of 11 burs designed by Dr. Fahrenkrug
- 4801931 Polishing Cups (RA)
- Pack of 5
 - Compatible with Contra Angle Handpiece (4804743)



4801447



4801931

Dental Burr Set (HP)


Tungsten carbide construction provides a long working life.

Includes:

- Round Burr – Size 18
- Flat Tapered Burr – Size 14 & 23
- Inverted Cone Burr – Size 12 & 21

Item No.

- 4802261 Dental Burr Set (HP)

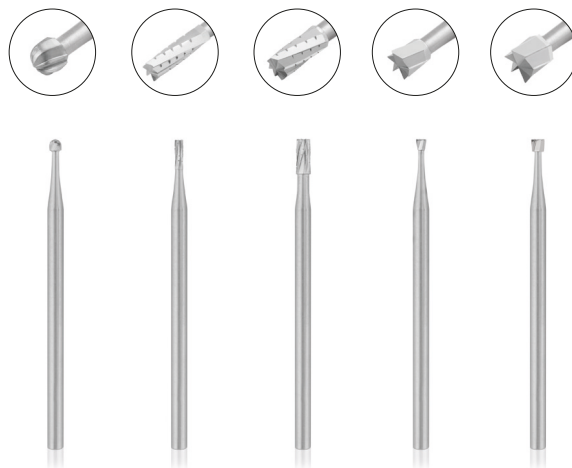



4800996



4804744



4804743

Low Speed Dental Handpieces

- Low speed dental handpieces are ideal to use for rodent dentals and dental polishing
- Slim design for easy handling
- Designed for air driven dental units
- Compatible with Profident (4808146) & Profident PLUS (4808148)

Item No.	
4800996	Low Speed Straight Handpiece with Air Motor <ul style="list-style-type: none">• 1:1, with internal water cooling• Up to 20,000 rpm• For operation of polishing attachments and cutting discs• Suitable for Profident Units
4804744	Low Speed Dental Handpiece – Straight (HP) <ul style="list-style-type: none">• Ideal for rodent dentals• Use with long straight shank burs, diamond discs and prophy heads• For polishing please order: Disposable Prophy heads HP (4808417)
4804743	Low Speed Dental Handpiece – Contra Angle (RA) <ul style="list-style-type: none">• For Dental Burr Set (4801447)• For polishing please order: Polishing Cups RA (4801931)

Eickemeyer®



*Innovations for
your everyday
success*

© Martin Carlssonc – stock.adobe.com

Disposable Polishing Heads (HP)

Perfect shiny teeth with a highly polished enamel surface.

Item No.

- | | |
|---------|---|
| 4808417 | Disposable Polishing Heads (HP) <ul style="list-style-type: none"> Single use Max. speed 3,000 U/min Compatible with Straight Angle (HP) Handpiece (4804744) Pack of 50 |
|---------|---|

**Polishing Cups (RA)**

Perfect shiny teeth with a highly polished enamel surface.

Item No.


- | | |
|---------|--|
| 4801931 | Polishing Cups (RA) <ul style="list-style-type: none"> Latch-type shank Compatible with Contra Angle Handpiece: 4804743 Pack of 5 |
|---------|--|


**Replacement bottle for Dental Units****Item No.**

- | | |
|---------|---|
| 4808151 | Replacement bottle for Dental Units <ul style="list-style-type: none"> Spare Water Bottle for Profident 2 litre |
|---------|---|

**Dental Drip Tray****Item No.**

- | | |
|---------|--|
| 4805230 | Dental Drip Tray <ul style="list-style-type: none"> Provides a comfortable position for the animal's head while helping to prevent water from soaking into a patient during scaling Removable grate for easy cleaning and disinfection Dimensions (in mm): L 390 x W 350 x H 45 |
|---------|--|

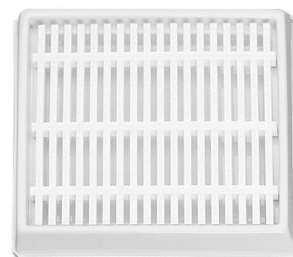


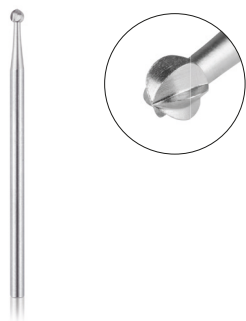

Dental Burr (HP)

Item No.

- 4809125 Dental Burr (HP)
- Size 23
 - Round

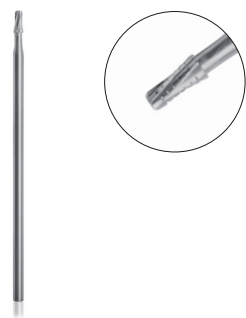


Dental Burr Set Rodent (HP)

Item No.

Length
mm

- | Item No. | | Length
mm |
|----------|-----------------------------|--------------|
| 4803173 | Dental Burr Set Rodent (HP) | 65 |
- Ideal for removing excess growth on molar teeth
 - Cross cut tapered fissure
 - Extra long shaft
 - Pack of 5

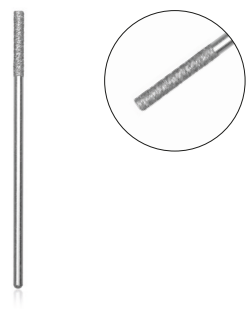

Dental Burr Rodent (HP)

Item No.

- 4802029 Dental Burr Rodent (HP)
- Ideal for removing excess growth on molar teeth
 - Diamond coated
 - For straight handpiece
 - Can be used with Dental Burr Guard Rodent (4803174)



Dental Burr Guard Rodent (HP)

The soft tissue protector guards the tongue and cheek during reshaping of occlusal surfaces.

Item No.

- 4803174 Dental Burr Guard Rodent (HP)

Diamond Disc Cutter SUPERFLEX

- Ideal for trimming rodent/rabbit incisor teeth
- Suitable for sectioning teeth
- Always wear eye protection (zz4809471)
- Diamond coated edge

Item No.		Diameter mm	Length mm
4802264	Diamond Disc Cutter SUPERFLEX (HP)	21	45
	<ul style="list-style-type: none"> • Use with Diamond Disk Safety Shield (4803175) 		
4805227	Diamond Disc Cutter SUPERFLEX (RA)	21	



4805227

Diamond Disc Safety Shield

- Protects surrounding soft tissue structures as well as the operator from spinning diamond disc
- Use with Diamond Disk (4802264)
- Always wear eye protection (zz4809471)

Item No.	
4803175	Diamond Disc Safety Shield

**Mouth Gags for Rodents**

An essential dental aid when examining the oral cavity of rabbits and other rodents.

Item No.		Length mm
4804244	Mouth Gag for Rodent	115
	<ul style="list-style-type: none"> • Suitable for chinchillas • Lightweight 	
4802096	Mouth Gag for Rodent	110
	<ul style="list-style-type: none"> • Suitable for medium to large breed rabbits • Robust 	



4804244




4802096

Mouth Gag for Rabbits

An essential dental aid when examining the oral cavity of rabbits, allows the veterinarian to have both hands free.

Item No.		Length mm
4803311	Mouth Gag for Rabbits	110
	<ul style="list-style-type: none"> • Cam ratchet mechanism (self-locking) • Robust 	



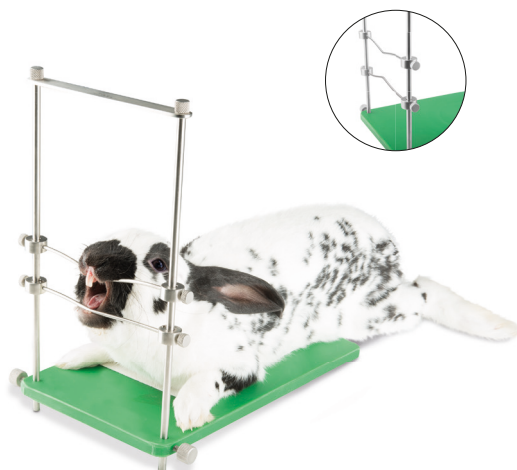

Rabbit/Rodent Tabletop Mouth Gag

The anaesthetised pet can be optimally fixed in order to gain the best possible access to the oral cavity. The worktable height and tilt are individually adjustable to achieve the optimal working position. The mouth opener is adjustable to fit any sized rabbit/rodent.

Item No.

- 4802194 Rabbit/Rodent Tabletop Mouth Gag
- Length (tabletop): 28 cm
 - Width (between the struts): 13.5 cm
 - Height: 29 cm

Rodent Cheek Dilators

Helps enhance examination of the oral cavity.

Item No.

Length mm

- | | | |
|---------|---|----|
| 4802093 | Rodent Cheek Dilator | 85 |
| | <ul style="list-style-type: none"> • Size: Medium • For medium to large rabbits • Blade dimension (in mm): L 21 x W 15 | |
| 4802094 | Rodent Cheek Dilator | 58 |
| | <ul style="list-style-type: none"> • Size: Small • For dwarf rabbits • Blade dimension (in mm): L 14 x W 11 | |
| 4802098 | Rodent Cheek Dilator | 58 |
| | <ul style="list-style-type: none"> • Size: Small • Long reach • Blade dimension (in mm): L 35 x W 10 | |
| 4803088 | Rodent Cheek Dilator | 60 |
| | <ul style="list-style-type: none"> • For chinchillas • Blade dimension (in mm): L 21 x W 9 tapering to L 21 x W 4 | |
| 4802058 | Rodent Cheek Dilator | 35 |
| | <ul style="list-style-type: none"> • Size: Extra Small • For small pets • Blade dimension (in mm): L 11 x W 9 | |



4802093



4802094



4802098



4803088



4802058



4801189



4801104

MINI FRIEDMANN Rongeurs

- Designed with horizontal grooved handles to enhance grip
- Double spring mechanism allows self-opening action


Item No.

Length mm


- | | | |
|---------|------------------------|-----|
| 4801189 | MINI FRIEDMANN Rongeur | 140 |
| | • Tip width: 1.3 mm | |
| 4801104 | MINI FRIEDMANN Rongeur | 120 |
| | • Tip width: 2.0 mm | |

Rodent Molar Cutter with Tungsten Carbide

Item No.		Length mm
4802100	Rodent Molar Cutter (TC) <ul style="list-style-type: none"> • Suitable for trimming long molars leaving teeth with a smooth surface • Double spring mechanism • Simple removal of hooks and tips of the molars 	170
4802817	Lewis Rasp – Fine Serration <ul style="list-style-type: none"> • Designed with back facing teeth 	180
4802818	Dental Rasp Lewis – Premium <ul style="list-style-type: none"> • Diamond coated 	180

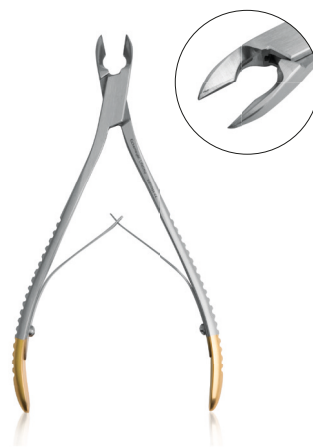

4802100

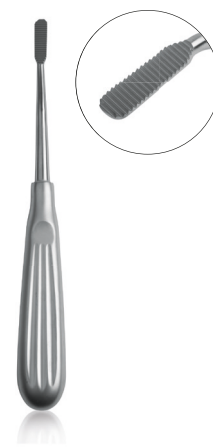

4802817

Rodent Molar Cutters

Helps to debulk spurs and hooks on molars while allowing visualisation of occlusal surfaces.

Item No.		Length mm
4802099	Rodent Molar Cutter <ul style="list-style-type: none"> • For chinchillas and rabbits • Spring mechanism 	180
4803310	Rodent Molar Cutter <ul style="list-style-type: none"> • For chinchillas • Angled tip • Spring mechanism 	120



4802099

*Quality equipment
for your practice*

Eickemeyer®



© Alex_Green – stock.adobe.com



4802093



4802094



4802098



4800234



4890236



4890237



4902092



4909557



4890240



4890241



4803173



4803174



4802246



4803175



4802095



4802096



4802097



4809780

Motorised Rodent Dental Kit

The most common dental disease in rodents is excessive tooth growth, which can affect both incisors and molars. Our motorised rodent dental kit has been designed by Richard Saunders (Veterinary Adviser to the Rabbit Welfare Association and Fund) and provides all the necessary instruments to help you treat commonly encountered problems.

Item No.		Length mm
4802237	Complete set, consisting of:	
4802093	Rodent Cheek Dilator <ul style="list-style-type: none"> • Size Medium • For medium to large rabbits • Blade dimension (in mm): L 21 x W 15 	85
4802094	Rodent Cheek Dilator <ul style="list-style-type: none"> • Size Small • For dwarf rabbits • Blade dimension (in mm): L 14 x W 11 	58
4802098	Rodent Cheek Dilator <ul style="list-style-type: none"> • Size: Small • Long reach • Blade dimension (in mm): L 35 x W 10 	58
4800234	Molar Extraction Forceps	160
4890236	CROSSLEY Incisor Luxator	195
4890237	CROSSLEY Molar Luxator	160
4802092	Mouth Gag for Rodent <ul style="list-style-type: none"> • Suitable for medium to large breed rabbits • Robust 	110
4809557	Peridontal Probe <ul style="list-style-type: none"> • Graduated tip 	155
4890240	Dental Cutters Rodent (HP) <ul style="list-style-type: none"> • Size 45 	50
4890241	Dental Cutters Rodent (HP) <ul style="list-style-type: none"> • Size 60 	50
4803173	Dental Burr Set Rodent (HP) <ul style="list-style-type: none"> • Ideal for removing excess growth on molar teeth • Cross cut tapered fissure • Extra long shaft 	65
4803174	Dental Burr Guard Rodent (HP)	
4802246	SUPERFLEX Diamond Disc Cutter (HP) <ul style="list-style-type: none"> • Ø 21 mm 	45
4803175	Diamond Disc Safety Shield	
4802095	Tongue Spatula <ul style="list-style-type: none"> • Straight • Flat 	175
4802096	HOLLMANN Tongue Spatula <ul style="list-style-type: none"> • Small • Concave 	190
4802097	HOLLMANN Tongue Spatula <ul style="list-style-type: none"> • Large • Concave 	190
4809780	Premium Instrument Box <ul style="list-style-type: none"> • Made of aluminium • Perforated bottom and lid • Features side closures • Not autoclavable 	



DENTANOMIC Dental Instruments

DENTANOMIC instruments are revolutionary dental instruments with blades that can be interchanged – making it easy for surgeons to switch to a different blade, or dispose of a blunt or damaged blade.

Features

- The light, ergonomically shaped handle is made of high-quality aluminium.
- The blades are especially hardened and have a sharp bevel to cut through the cervical periodontal fibres with minimal effort and widen the periodontal gap.

Item No.		Width mm
4802541	DENTANOMIC Ergonomic Handle • Blue	
4802542	DENTANOMIC Ergonomic Handle • Green	
4802387	DENTANOMIC Ergonomic Handle • Pink	
4802551	DENTANOMIC Ergonomic Handle • Longer shaft • Orange	
4802463	DENTANOMIC Ergonomic Handle • Longer shaft • Purple	
4802554	DENTANOMIC Winged Elevator Tip	2
4802056	DENTANOMIC Winged Elevator Tip	3
4802555	DENTANOMIC Winged Elevator Tip	4
4802556	DENTANOMIC Winged Elevator Tip	5
4802553	DENTANOMIC Luxation Tip	1.5
4802557	DENTANOMIC Luxation Tip	2
4802558	DENTANOMIC Luxation Tip	3
4802559	DENTANOMIC Luxation Tip	4
4802065	DENTANOMIC Luxation Tip	5
4802560	Storage and Sterilisation Rack • For 9 DENTANOMIC blades • With lid • Dimensions (in cm): L 19 x W 8 x H 2	
4805652	DENTANOMIC Honing Pack • 1 DMT DiaSharp (extra-fine) & 1 DMT DiaFold sharpener (extra-fine) • Diamond Hone	

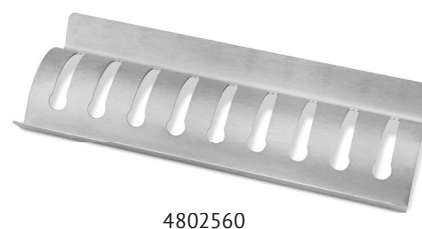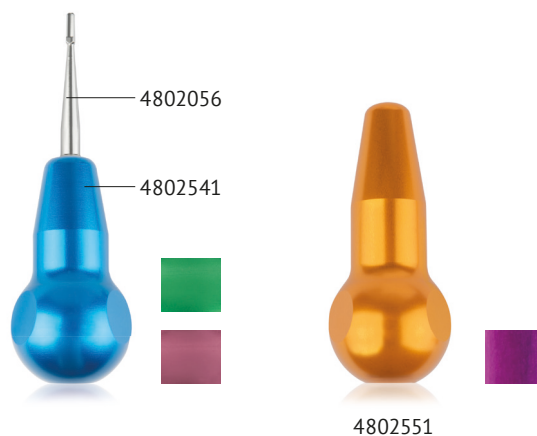

Winged Elevators

Designed with a sharp blade and lateral wings, which improves the contact with conical roots.

Item No.		Width mm	Length mm
4802251	Winged Elevator	2	145
4800644	Winged Elevator	3	145
4800645	Winged Elevator	5	145





4802086

BEIN Dental Elevators

The spoon shaped configuration of these elevators helps guide them into the periodontal pocket, which gives additional leverage and reduces gingival damage.

Item No.		Width mm	Length mm
4802086	BEIN Dental Elevator	2	135
4802087	BEIN Dental Elevator	3	135
4802088	BEIN Dental Elevator	4	135




4808460

BEIN Dental Elevators Modified

The spoon shape configuration of these elevators helps guide them into the periodontal pocket, which provides additional leverage and reduces gingival damage.

Item No.		Width mm	Length mm
4808458	BEIN Dental Elevator Modified	2	135
4808459	BEIN Dental Elevator Modified	3	135
4808460	BEIN Dental Elevator Modified	4	135

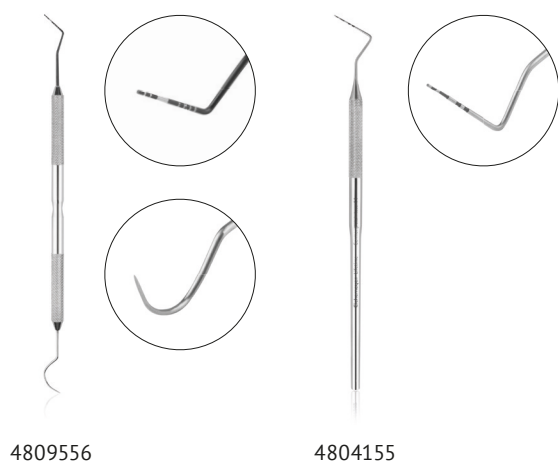

4809556

4804155

Periodontal Probes

Designed with a blunt ended tip for subgingival exploration.

Item No.		Length mm
4809556	Periodontal Probe • Dual-ended • Graduated tip	170
4804155	Periodontal Probe • Graduated tip	160
4801416	Periodontal Probe/Explorer • Size 9	155
4804156	Periodontal Probe/Explorer • Size 17	170

Dental Instruments for Small Animals

Item No.		Length mm
4804152	JAQUETTE Supragingival Scaler	165
4804153	Subgingival Curette <ul style="list-style-type: none"> Designed with a 70° blade to the shaft 	165
4801438	Dental Mirror <ul style="list-style-type: none"> Width: 22 mm 	165



4804152




4804153



4801438

Subgingival Curettes/Dental Periosteal Elevators

Item No.		Length mm
4808780	Subgingival Curette Gracey <ul style="list-style-type: none"> Size 3/4" Suitable for incisor teeth 	170
4808786	Subgingival Curette Gracey <ul style="list-style-type: none"> Size 9/10" Suitable for buccal and lingual aspects of molar teeth 	170
4801202	WILLIGER Periosteal Elevator <ul style="list-style-type: none"> Width: 6.5 mm 	160

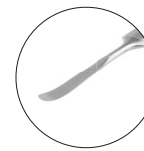

4808780



4808786



4801202

FAHRENKRUG Dental Elevator


Designed to help extract deciduous teeth.

Item No.		Width mm	Length mm
4804145	FAHRENKRUG Dental Elevator <ul style="list-style-type: none"> For dogs 	2	200
4801314	FAHRENKRUG Dental Elevator <ul style="list-style-type: none"> For cats and smaller dogs 	3	190



4804145



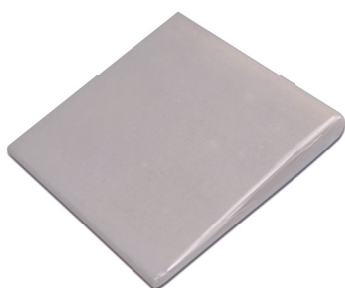

Arkansas Whetstone

- Suitable for sharpening dental elevators
- Use oil as a lubricant to prevent steel particles from clogging the stone pores

Item No.		Width cm	Length cm
4800595	Arkansas Whetstone	5.1	5.1



Dental Extraction Forceps

There are various extraction forceps configurations to “fit” different tooth shapes and should be used alongside dental elevators to help remove teeth.

Item No.		Length mm	Form
4805144	Dental Extraction Forceps • Diamond coated	110	
4802089	Dental Extraction Forceps • Premolars	110	straight
4802090	Dental Extraction Forceps • Roots	165	straight
4805144	Dental Extraction Forceps • Premolar, incisors, roots	165	straight



WITZEL Root Extraction Forceps

Designed to extract root fragments from the maxillary arcade.

Item No.		Length mm
4801414	WITZEL Root Extraction Forceps	150

STIEGLITZ Root Extraction Forceps


Designed with narrow tips to help grasp fractured tooth roots in cats and small breed dogs.


Item No.		Length mm
4808250	STIEGLITZ Root Extraction Forceps	130

**Dental Hand Scalers**

Helpful to remove dental tartar.

Item No.	
4804146	Dental Hand Scaler • Hook shaped
4804147	Dental Hand Scaler • Enamel hatchet
4804148	Universal Hand Scaler • Left
4804149	Universal Hand Scaler • Right
4804150	Dental Hand Scaler • Hatchet
4804151	Dental Hand Scaler • Hatchet • Robust

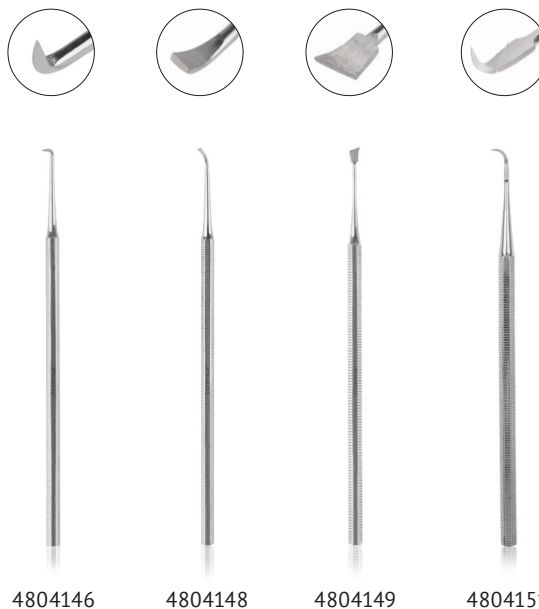

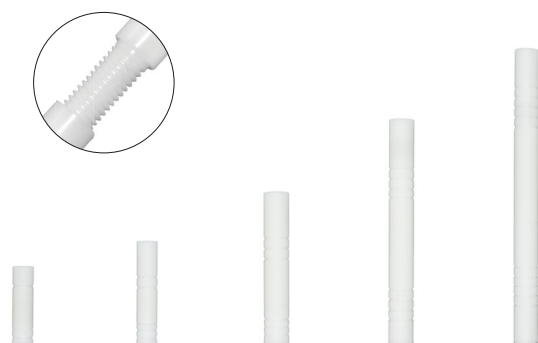

EICKEMEYER® Mouth Gag Set according to KOJIMA

With the EICKEMEYER® mouth gag set according to Kojima dental procedures and dental corrections in cats and dogs can now be carried out even more secure. These novel mouth gags consist of two cylinders that can be adjusted via a thread in their length. With this feature the mouth gags can be easily adapted to the mouth of the patient. By the defined adjustment an excessive tension is no longer possible, therefore irreparable damages and injuries can be avoided.

- Material: Hard plastic
- Easy to clean and disinfect
- Radiolucent
- Set of 5 sizes

Item No.

- 4801146 Complete set, consisting of:
- Size XS (3 – 4 cm)
 - Size S (4 – 6 cm)
 - Size M (6 – 8 cm)
 - Size L (9 – 11.5 cm)
 - Size XL (12 – 15 cm)

Wood Dental Mouth Gags

Being made of wood makes these mouth gags exceptionally lightweight and allows secure access to the mouth.

Item No.		Length mm
4801567	Wood Dental Mouth Gag • Small	70
4801568	Wood Dental Mouth Gag • Large	90



4801567

4801568

Dental Models

- Transparent hinged jaw helps visualise the roots of all teeth (incisors, canine, premolars, molars)
- Can be separated for closer examination
- Ideal for client education

Item No.

- 4804479 Dental Model Cat
- 4803362 Dental Model Dog
- 4802233 Dental Model Rodent



4803362




4802233

Canine and Feline Dental Instrument Kit

Our dental kit provides all the necessary instruments to help you tackle the most commonly encountered dental conditions. It includes the new EICKEMEYER® Mouth Gag Set according to KOJIMA, so dental procedures and dental corrections in cats and dogs can now be carried out even more securely.

Item No.		Length mm	Form
4801512	Complete set, consisting of:		
4808248	Apical Elevator	145	straight
4802251	Winged Elevator • Width: 2 mm	145	
4808458	BEIN Dental Elevator Modified • Width: 2 mm	135	straight
4808459	BEIN Dental Elevator Modified • Width: 3 mm	135	straight
4808460	BEIN Dental Elevator Modified • Width: 4 mm	135	straight
4804145	BEIN Dental Elevator • For dogs • Width: 2 mm	200	
4801314	BEIN Dental Elevator • For cats and smaller dogs • Width: 3 mm	190	
4809556	Periodontal Probe • Dual-ended • Graduated tip	170	
4808780	Subgingival Curette Gracey • Size 3/4" • Suitable for incisor teeth	170	
4808786	Subgingival Curette Gracey • Size 9/10" • Suitable for buccal and lingual aspects of molar teeth	170	
4808250	STIEGLITZ Root Extraction Forceps	130	
4802089	Dental Extraction Forceps • Premolars	110	straight
4801146	EICKEMEYER® Mouth Gag according to KOJIMA Complete set, consisting of: • Size XS (3 – 4 cm) • Size S (4 – 6 cm) • Size M (6 – 8 cm) • Size L (9 – 11.5 cm) • Size XL (12 – 15 cm)		
4809780	Premium Instrument Box • Made of aluminium • Perforated bottom and lid • Features side closures • Not autoclavable • Dimensions (in cm): L 30 x W 14 x H 7		

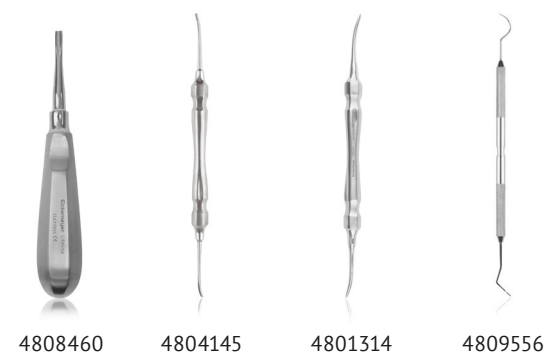



4808285



4809645



4805841



4802155

DENTO-DORM® Positioning Accessories

- DENTO-DORM® products are developed, produced and clinically tested according to strict guidelines in Germany
- Warm insulating positioning aid constructed of sturdy hard foam
- Suitable for dental procedures, X-Rays and other surgical procedures for dogs, cats and small pets

Item No.

- | | |
|---------|--|
| 4808285 | DENTO-DORM® Mattress <ul style="list-style-type: none"> • Includes: Four positioning wedges • Dimensions (in cm): L 75 x W 30 x H 8 |
| 4809645 | DENTO-DORM® Head Positioning Aid <ul style="list-style-type: none"> • Two per set |

Low Lift Scissor Table II

- Mobile Scissor lift examination and surgical table
- Foot pedal operated battery electric lift to get the correct working height and protect your back
- Adjustable from 28 cm to 110 cm high
- Smooth stainless steel flat top with anti-drip edges and drain hole and wooden base for reduced noise
- Fixed tie down cleats
- Table top 60 cm x 130 cm
- Includes mobile locking, swivelling caster kit (4802154), battery charger (4808433), tie down kit (4802303), battery (4800096), lift motor (4806823), foot control (4808434)

Item No.

- | | |
|---------|---|
| 4805841 | Low-Lift Scissor Table II |
| 4802155 | Dental Tray for Scissor Table II <ul style="list-style-type: none"> • Easy to clean and disinfect stainless steel • Pluggable drain • Dimensions (in cm): W 48.5 x H 5.5 x L 35.5 |

DR. MACH LED 130 F

DR. MACH from Germany has been setting the standards in medical illumination technology for decades. DR. MACH's most recent version of our most popular LED light is the LED 130 F. LED 130 F offers 70,000 Lux light intensity with a focusable beam from the conveniently placed handle. The biggest advantages of LED technology is the guaranteed minimum life-span of 40,000 hours, energy efficiency and virtually no heat operation.

- 70,000 Lux at 1 m distance
- Colour rendering index Ra \geq 95
- Dimmable; 100 % to 50 %
- Light head diameter 33 cm
- Light field 15 – 21 cm
- Number of LEDs: 19
- Power consumption: 28 W
- Include 110 V/60 hZ electrical plug (4809019)

Item No.

4801291	DR. MACH LED 130 F <ul style="list-style-type: none"> • Mobile floor stand 4' / 120 cm high with SWING arm
4803659	DR. MACH LED 130 F <ul style="list-style-type: none"> • Wall mount
4801292	DR. MACH LED 130 F <ul style="list-style-type: none"> • Ceiling mount; for ceilings up to 9' / 3 m
4809007	Sterilizable light handle conversion kit <ul style="list-style-type: none"> • Replaces existing handle with removable autoclavable handle



4801291



4803659



4801292

CLEANING & CARE

Dental Handpieces

With regular care and cleaning after each procedure, you can maintain the full functionality of your handpieces.

1. Preparation

After each procedure, remove the handpiece from the coupling and then remove the drill (ideally with a pair of tweezers).

2. External cleaning of the handpieces/turbines

The handpiece should be cleaned with warm, running water (e. g. with a toothbrush) so that abrasive dust or saliva is washed away. Never clean the handpieces in an ultrasonic bath!

3. Disinfection of handpieces/turbines

The handpieces can be treated externally in the spray water/spray air channels with disinfectants. Never immerse handpieces in disinfectant solution!

4. Care of handpieces/turbines

Handpieces should be sprayed with care oil after each clean. Spray the care oil into the water and spray channels of the handpieces for 1 – 2 seconds. If the liquid emerging at the head is still mixed with dirt, the process should be repeated.



Injecting maintenance oil into the spray channel (maintenance oil item no. 4801326 and spray channel adapter item no. 4801328)



Injecting care oil into handpieces with a planar receiving shaft (care oil item no. 4801326 and handpiece adapter item no. 4800420)

5. Care of collet chucks

The free movement of the push buttons is very important for the fixed clamping function of the drill. The clamping forceps must be sprayed with the maintenance oil at least once a week from the drill side. Only use diamond grinders and drills from well-known manufacturers. If the shaft material is too soft, scrapings may occur if the drill is twisted in the clamping jaw. This results in fine chips, which can become caught in the clamping system. Thereafter, the drill may no longer be removed and the tensioning system thereby damaged.

6. Care

It is important to keep the all instrument joints well lubricated with instrument oil (e. g. no. 4803909) to avoid any friction or corrosion. Other care checks include smoothness of movement. We offer a repair service for instruments.

OIL REGULARLY FOR A LONG LIFESPAN

CONTENT**A**

Arkansas Whetstone..... 16

C

Cleaning and Care for Dental Handpieces 22

D

Dental Burr Guard.....8

Dental Burrs (FG).....4

Dental Burrs (HP)4, 5, 8

Dental Burrs (RA).....5

Dental Drip Tray7

Dental Hand Scalers.....17

Dental Handpieces4, 6

Dental Instruments for Small Animals 13, 15

Dental Models 18

Dental Unit Air Compressor4

Dental Units..... 2, 3

Dentistry Kits for Dogs and Cats..... 19

Dentistry Kits for Rodents..... 12

Diamond Disc Cutter SUPERFLEX.....9

Diamond Disc Safety Shield.....9

E

Elevators 13–15

Examination and Theatre Lights..... 21

Extraction Forceps 16, 17

M

Mouth Gags for Rodents9, 10

Mouth Gags for Small Animals 18

P

Periodontal Probes.....14

Polishing Cups (RA) 7

Polishing Heads (HP) 7

Positioning Accessories..... 20

R

Replacement bottle for dental units 7

Rodent Cheek Dilators.....10

Rodent Molar Cutters.....11

Rongeurs.....10

S

Scissor Table.....20

Subgingival Curettes15

EICKEMEYER Inc.
617 Douro Street, Suite #205
Stratford, Ont. Canada
N5A 0B5
T +1 519 273 5558
F +1 519 271 7114
info@eickemeyervet.ca
www.eickemeyer.com